

CPAP MACHINE SUPERSTAR



ELK-200B



ELK-200A



ELK-200C



**ELK-200D
Bubble**

ELK-200B CPAP System

Application

Continuous Positive Airway Pressure (CPAP) nurtures and supports infant breathing by providing respiratory support throughout the respiratory cycle. CPAP System maintains the infant's functional residual capacity by helping to prevent airway closure. ELK-200B CPAP System features an ergonomic and uncluttered touch-screen interface, enabling the clinician to quickly and easy access the current patient status. With the advantages of compact device, small footprint, saving energy and integrated high-performance in operation. This makes it an ideal companion for your smallest patients in the Neonatal Intensive Care Units (NICU), emergency ward, delivery room and other departments.

TRUST POINT

- 12.1 Inch TFT touch LCD screen displays the ventilation parameters, alarm information and oscillogram.
- Better work of breathing: Integration breathing circuit design, maximize lung protection and breathing support.
- Increase clinician's confidence in Neo-CPAP: Multiple parameters monitoring interface, makes every parameter clear, let users know the patient conditions in all aspects.
- Real time pressure-time oscillogram and high precision O₂ concentration detection function included, reduces length of stay in hospital
- Built-in air and oxygen blender, ensures stable oxygen concentration.
- Stable and low noise air compressor creates quiet work environment for doctors, can be used for 2-3 Neo-CPAP Systems together.



ELK-200B CPAP System



Intuitive operation

12.1 Inch LCD touch screen displays measured parameters clearly and Intuitively

Minimize noise effect

Reducing noise to 50dB during using, to keep the ward quiet



Minimizes airway drying and inflammation with intelligent humidification

Adjustable temperature and humidity settings, perfectly adapt the humidification settings to the special needs of little patients and environmental conditions. Total 3 types for choice: Thermometer, Digital, Heating wire



Multilevel filtration-purification system

Reduce in maximum the dust effect, increase of service life

Technologically advanced safety valve

Improve lung volume recruitment and maintain infant energy reserves, without the complications associated with endotracheal intubation



Smart backup battery

power source goes off, backup power source starts to work automatically, and last 2-3 hours



Mobility for intrahospital transport

Easy to move with 4 casters, easy to stop with 2 brakes in front casters



TECHNICAL SPECIFICATION

Ventilation Mode

CPAP, NCPAP, NIPPV, HFNC, MANU

Adjusting Parameter

Oxygen concentration	21 ~ 100%
CPAP	0 cmH ₂ O ~ 15 cmH ₂ O (NCPAP)
Continuous flow	1 L/min ~ 20 L/min (HFNC) 2 L/min ~ 20 L/min (CPAP)
Manual ventilation continuous time	1 ~ 15s
Manual ventilation pressure	3 ~ 15 cmH ₂ O (NIPPV, NCPAP)
Manual ventilation flow	3 ~ 25 L/min (HFNC, CPAP)
Fast continuous O ₂ supply	30 ~ 120s
Fast O ₂ supply concentration	23% ~ 100%
IPAP	3 cmH ₂ O ~ 20 cmH ₂ O (NIPPV)
EPAP	1 cmH ₂ O ~ 13 cmH ₂ O (NIPPV)
Breath Frequency	1 /min ~ 120 /min (NIPPV)
Adjusting Pressure Limit	1 cmH ₂ O ~ 60 cmH ₂ O Lower: 0 cmH ₂ O ~ 15 cmH ₂ O

Monitoring Parameter

Oxygen concentration	15 ~ 100%
IPAP/EPAP	0 cmH ₂ O ~ 60 cmH ₂ O (NIPPV)
Airway pressure	0 cmH ₂ O ~ 60 cmH ₂ O (NCPAP, HFNC)
Continuous flow	0 L/min ~ 20 L/min
Breath frequency	0 /min ~ 150 /min
Inhalation time	0.1 s ~ 20 s
ETCO ₂ (Optional)	0 ~ 152 mmHg

Oscillogram

P-T(pressure-time)

Oxygen gas flow column

Air gas flow column

Air oxygen mixed gas flow column

Alarm and Protection

Low Air pressure alarm	≤ 0.2 MPa
Low O ₂ pressure alarm	≤ 0.2 MPa
High oxygen concentration alarm	19% ~ 100%
Low oxygen concentration alarm	18% ~ 99%
High Airway pressure alarm	1 cmH ₂ O ~ 60 cmH ₂ O
Low Airway pressure alarm	0 cmH ₂ O ~ 15 cmH ₂ O
Internal battery backup low voltage alarm	<10.2 V ± 0.3 V
Pressure imbalance alarm	The Air, O ₂ input pressure > 200kPa, And Air, O ₂ input pressure difference ≥ 100 kPa.
Maximum limit pressure alarm	≤ 60 cmH ₂ O
Air compressor alarm	≥ 100 °C ± 10 °C
AC power failure alarm	
Fan failure alarm	
Silent alarm	

Working Condition

Voltage	100 ~ 240V
Frequency	50/60 HZ
Input power	900 VA (with air compressor) 60 VA (without air compressor)
Gas source	O ₂ , AIR (medical grade)
Pressure	350 kPa ± 50 kPa
Ambient temperature range	+ 5°C ~ + 40°C
Relative humidity range	≤ 80%
Atmospheric pressure range	860hPa ~ 1060hPa

Wooden Case Packing Size

Main components : L 480 * W 495 * H 580 mm

G.W. : 22 KG N.W. : 10 KG

Air compressor : L 683 * W 687 * H 140 mm

G.W. : 136 KG N.W. : 54 KG

ELK-200A CPAP System

Application

Continuous Positive Airway Pressure (CPAP) nurtures and supports infant breathing by providing respiratory support throughout the respiratory cycle. CPAP System maintains the infant's functional residual capacity by helping to prevent airway closure.

The ELK-200A CPAP System has the advantages of compact device, small footprint, saving energy and integrated high-performance in operation. This makes it an ideal companion for your smallest patients in the Neonatal Intensive Care Units (NICU), emergency ward, delivery room and other departments.



TRUST POINT

- Easy-to-understand monitoring: The configurable screen displayIntuitionistic simplicity and specified approach diverse monitoring data.
- Mobility for intrahospital transport: Easy to move with 4 wheels, easy to stop with two brakes in front casters.
- Minimizes airway drying and inflammation with Intelligent humidification: Adjustable temperature and humidity settings, perfectly adapted the humidification settings to the special needs of small patients and environmental conditions.
- Minimize noise effects: Reducing noise to 50dB during using, to keep the ward quiet.
- Designed and manufactured by Superstar Medical with over 25 years experience in this area, Over 2,000 units installed in the world.

FEATURE

- 7" TFT LCD screen displays the ventilation parameters, alarm information and oscillograms.
- Intergrated breathing circuit design, ensure easy operating and keep hygiene.
- Multiple parameters monitoring interface ensures every parameter is clear, let operator know the patient conditions in all aspects.
- Real time pressure-time oscillogram and high precision oxygen concentration detection function included.
- Built-in air and oxygen blender, ensure stable oxygen concentration.
- Stable and low noise air compressor creates quiet work environment for doctors and patients, able to support 2-3 CPAP machines at the same time.

SAFETY


- Three levels alarm system, visual and audible alarm information.
- Different types of alarm, reminder and protection functions are available.
- Advanced power management control technology.
- Built-in backup battery provides the emergency power supply to the unit.



7" TFT LCD screen



Shuttle button

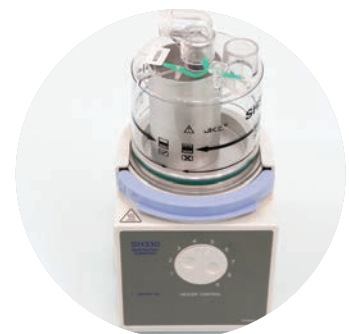

Humidifier

TECHNICAL SPECIFICATION

Ventilator Parameter Range

Oxygen concentration	Adjusting range: 21% - 100% Monitoring range: 15% - 100%
CPAP	Adjusting range: 0cmH ₂ O - 15cmH ₂ O Monitoring range: 0cmH ₂ O - 60cmH ₂ O
Continuous flow	Adjusting range: 2L/min - 20L/min Monitoring range: 0L/min - 20L/min
Manual ventilation	Output flow: 20L/min ± 5L/min
Air compressor	Continuous output flow: ≥ 25L/min Peak output flow: ≥ 120L/min
Oscillogram	P-T (pressure - time)

Working condition

Gas source	O ₂ , Air
Pressure	280kPa - 600kPa
Voltage	100 - 240V
Power frequency	50/60Hz

Alarm and Protection

AC power failure alarm	Power failure or no connection
Internal backup battery low voltage alarm	< 11.3 ± 0.3V
High oxygen concentration alarm	19% - 100%
Low oxygen concentration alarm	18% - 99%
High airway pressure alarm	1.5kPa - 6kPa
Low air input pressure alarm	< 200kPa
Low oxygen input pressure alarm	< 200kPa
Air and oxygen blender alarm	When any gas source lack When differential pressure > 0.1mPa
Compression pump alarm	Excess temperature

Size

Main engine	L 480*W 670*H 450mm
G.W.	25KG
N.W.	12KG
Air compressor	L 683*W 687*H 1140mm
G.W.	94KG
N.W.	63KG

ELK-200C CPAP System

Application

Continuous Positive Airway Pressure (CPAP) nurtures and supports infant breathing by providing respiratory support throughout the respiratory cycle. CPAP System maintains the infant's functional residual capacity by helping to prevent airway closure.

The ELK-200C CPAP System has the advantages of compact device, small footprint, saving energy and integrated high-performance in operation. This makes it an ideal companion for your smallest patients in the Neonatal Intensive Care Units (NICU), emergency ward, delivery room and other departments.



TRUST POINT

- Easy-to-understand monitoring: The configurable screen display, Intuitionistic simplicity and specified approach diverse monitoring data.
- Mobility for intrahospital transport: Easy to move with 4 wheels, easy to stop with two brakes in front casters.
- Minimizes airway drying and inflammation with Intelligent humidification: Adjustable temperature and humidity settings, perfectly adapted the humidification settings to the special needs of small patients and environmental conditions.
- Minimize noise effects: Reducing noise to 50dB during using, to keep the ward quiet.
- Designed and manufactured by Superstar Medical with over 25 years experience in this area, Over 2,000 units installed in the world.



3.5" TFT LCD screen



Shuttle button



Humidifier

Technical specification

Oxygen concentration	Adjusting range: 21% - 100% Monitoring range: 15% - 100%
CPAP	Adjusting range: 0cmH ₂ O - 15cmH ₂ O Monitoring range: 0cmH ₂ O - 60cmH ₂ O
Continuous flow	Adjusting range: 2L/min - 20L/min Monitoring range: 2L/min - 20L/min
Gas source	O ₂ , Air
Pressure	280kPa - 600kPa

Voltage	100 - 240V
Power frequency	50/60Hz
Main engine	L 480*W 670*H 450mm
G.W.	25KG
N.W.	10KG
Air compressor	L 683*W 687*H 1140mm
G.W.	94KG
N.W.	63KG

ELK-200D Bubble CPAP System

Application

Continuous Positive Airway Pressure (CPAP) nurtures and supports infant breathing by providing respiratory support throughout the respiratory cycle.

CPAP System maintains the infant's functional residual capacity by helping to prevent airway closure.

The ELK-200D CPAP System has the advantages of compact device, small footprint, saving energy and integrated high-performance in operation. This makes it an ideal companion for your smallest patients in the Neonatal Intensive Care Units (NICU), emergency ward, delivery room and other departments.



TRUST POINT

- Mobility for intrahospital transport : easy to move with 5 wheels, easy to stop with three brakes in front casters.
- Flexible configurations to suit your needs.
- International standard and advanced technology suitable for all baby patient around the world.
- Designed and manufactured by Superstar Medical with over 25 years experience in this area, Over 2,000 units installed in the world.



Flowmeter



Humidifier



Pressure gauge

Feature

- Integrated breathing circuit design, ensure easy operating and keep hygiene.
- Built-in air and oxygen blender, ensure stable oxygen concentration.
- Reliable CPAP control valve and pressure monitoring system improve CPAP adjustment precision.

Safety

- Two casters with brake, ensure stable operating.
- Pressure gauge and flow meter ensure accurate adjusting of parameters.

TECHNICAL SPECIFICATION

Oxygen concentration	Adjusting range: 21% - 100%
CPAP	Adjusting range: 0cmH ₂ O - 10cmH ₂ O
Continuous flow	Adjusting range: 1.5 L/min - 15 L/min
Gas source	O ₂ , Air
Pressure	280kPa - 600kPa

Wooden case packing size (main component)	L 700*W 710*H 380mm
---	---------------------

G.W.	30KG
------	------

N.W.	17.7KG
------	--------

Wooden case packing size (Trolley)	L 160*W 195*H 1100mm
------------------------------------	----------------------

G.W.	8KG
------	-----

N.W.	2KG
------	-----

Total CBM	0.7m ³
-----------	-------------------