

Dräger Bubble CPAP Therapy Total Solution Consumables and Accessories

The Seattle-Positive Airway Pressure (PAP) system is an innovative solution which uses the proven advantages of Bubble CPAP therapy - such as oscillatory effects like high frequency ventilation - and combines it into a unique design enabling easy access to affordable healthcare. The total Dräger solution for Bubble CPAP is made up of durable and high quality products providing an affordable solution with increased therapeutic value to help infants with respiratory distress breathe easier.



Infusion support

Oxymixer Air and O₂ blender


Hinge Arm

Air and O₂ hoses

Bubble CPAP circuit -
Seattle PAP

Humidifier

Trolley

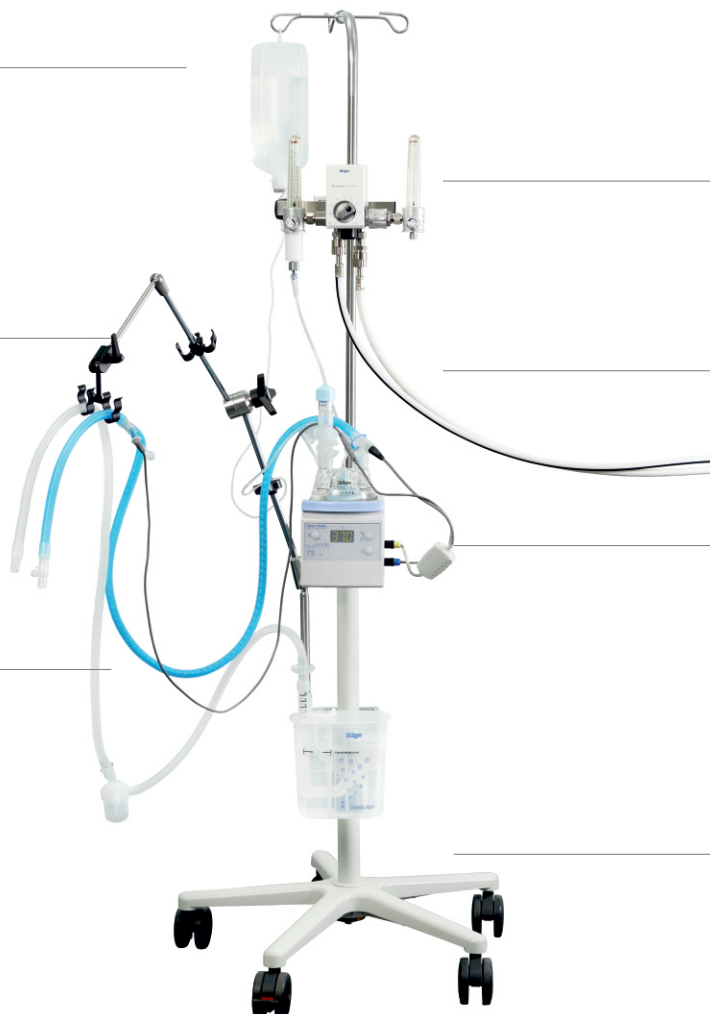

Benefits

Answering the Challenge

Worldwide, almost 1 million infants die each year from respiratory distress.¹ Respiratory muscle exhaustion is a major cause of failure in conventional n-B-CPAP therapy. These failures, which occur in half of babies receiving this therapy, often necessitate endotracheal intubation and mechanical ventilation that carry the risk of tissue damage and inflammation.

The Seattle-PAP from Dräger may be a solution for such complications. The patented Seattle-PAP Bubble CPAP System is designed to create continuous positive airway pressure and pressure oscillations. The 135° angle in the submersed tubing increases the amplitude of oscillations, thus providing effects similar to high frequency oscillatory ventilation for improving gas exchange and offer more efficient respiratory support compared to other methods of noninvasive respiratory support. These oscillations are thought to lessen the work required from respiratory muscles, which would make it easier for babies to breathe.

Regrettably, most of these casualties occur in resource-limited countries due to a lack of access of respiratory support devices commonly used in the developing world. We at Dräger take responsibility for those who can't afford high-end technology. The Seattle-PAP has the potential to make effective respiratory support accessible to newborns anywhere in the world.

Easy to set-up, easy to use

Seattle-PAP is safe, effective and easy to use. It was designed to:

- Promote infant safety with a uniquely designed pressure manifold with pressure relief valve.
- Safety tube lock to avoid unintended changes of the PEEP level.
- Provide comfort and optimum CPAP support to the infant due to disposable BabyFlow and BabyFlow plus interfaces complementing the system.
- Include funnel with the filling tube which minimizes spillages.
- The swivel connector to the expiratory limb prevents unnecessary disconnections and pressure loss.
- The water trap in the expiratory limb prevents unwanted pressure peaks due to condensation in the expiratory limb.

It requires only oxygen and compressed air, both to be administered by an O2-blender such as the Oxymixer from Dräger and a humidification source for the circuit.

Simple yet versatile

The Seattle PAP has a functional design with a complete focus on the infant's clinical needs. It is both easy to understand and use. You can:

- Easily adjust the CPAP from 4.5 to 10 cm H₂O.
- Quickly connect the Bubble CPAP System to your choice of BabyFlow CPAP interface solutions.
- In combination with BabyFlow plus you can treat patients with the needed masks and prongs.

Benefits

Safety first

Inadvertent changes of the position of the tube in the water, thus changing the CPAP level can cause serious harm to your patient. With the Seattle PAP you select your desired settings by locking the tube in the set position. This avoids accidental changes of CPAP levels and is a unique safety feature. Safety and simplicity with no need to worry.


Economical Technology

New innovations and technology should also be affordable for developing countries. Due to the simplicity of the design, the Seattle-PAP is exactly that – affordable respiratory care. This means that effective and life saving therapy is now available to many parts of the world that did not have access before.

Four versions for maximum flexibility with our Blenders

The Oxymixer is available in four different versions: a high-flow version and a low-flow version, each of which can be supplied with or without oxygen monitoring. The Oxymixer comes complete with flow meters which control the flow of dosage. All versions allow a high degree of flow accuracy. The low-flow version Oxymixer with bleed even allows a highly accurate flow with rates of under 3 l/min.

Details

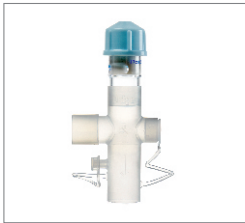

Seattle PAP and Seattle PAP Plus each pack contains:

Inspiratory heated neonatal circuit, expiratory watertrap

Bubble CPAP chamber

Humidifier chamber

Safety pop off valve

Oxygen hose

Blender Oxymixer



Oxymixer low flow without monitor



Oxymixer low flow with silencer



Oxymixer low flow with monitor

Accessories

D-10096-2018



BabyFlow plus

The superior non-invasive respiratory support system, designed to maximize performance and patient comfort. Adjustable angle to ensure best possible fit of masks and prongs.

D-10096-2018



Prongs

The nasal prongs are available in 11 sizes and two different widths. In combination with the BabyFlow plus, the prongs deliver an effective seal and reduce pressure on the infant's septum, therefore, minimizing skin breakdown.

D-10097-2018



Masks

The nasal masks are available in six different sizes. Sizes S, M and L are well known and accepted silicone masks for smaller patients up to approximately 2 kg (4.41 lbs). The new masks for patients from 2 kg (4.41 lbs) up to 5 kg (11 lbs) are designed for maximum therapeutic effectiveness and patient comfort using 3D infant facial scan technology. The flexible bellows allow the mask to find a natural position over the infant's nose, which helps minimize leaks and maintain a consistent nCPAP level.

D-10098-2018



Headgear

The headgear helps secure the BabyFlow plus quickly and easily while minimizing pressure points. It can be gently wrapped around the infant's head and secured with adjustable closures. It is available in six sizes to fit head circumferences from 17 cm to 42 cm (6.69 inches to 16.54 inches).

D-10099-2018



Microfiber Caps

The flexible and breathable microfiber caps are available in seven sizes. The disposable caps include side straps and top straps to help secure and stabilize the BabyFlow plus.

D-100100-2018



Self Connector

Enables the connection to the majority of neonatal breathing circuits.

Ordering Information

Seattle PAP		MP02250	
Seattle PAP plus		MP02260	
BabyFlow® plus disposable			
S	20 pieces		MP03700
M	20 pieces		MP03701
L	20 pieces		MP03702
Microfiber Caps			
XS	white	5 pieces	84 18 533
S	yellow seam	5 pieces	84 18 534
M	red seam	5 pieces	84 18 535
L	green seam	5 pieces	84 18 536
XL	orange seam	5 pieces	84 18 537
XXL	blue seam	5 pieces	84 18 538
XXL plus	white/black seam	5 pieces	84 18 539
Headgear			
XXS	5 pieces		MP03880
XS	5 pieces		MP03881
S	5 pieces		MP03882
M	5 pieces		MP03883
L	5 pieces		MP03884
XL	5 pieces		MP03885
Neo mask			
S	10 pieces		84 18 491
M	10 pieces		84 18 490
L	10 pieces		84 18 619
4	10 pieces		MP01494
5	10 pieces		MP01495
6	10 pieces		MP01496
Standard prongs			
Prongs XS	10 pieces		84 18 415
Prongs S	10 pieces		84 18 605
Prongs M	10 pieces		84 18 416
Prongs L	10 pieces		84 18 531
Prongs XL	10 pieces		84 18 417
Prongs XXL	10 pieces		84 18 617
Wide prongs			
Prongs 3	10 pieces		84 18 603
Prongs 4	10 pieces		84 18 607
Prongs 5	10 pieces		84 18 609
Prongs 6	10 pieces		84 18 611
Soft connector N	20 pieces		MP03826
Demo pack BabyFlow® plus, size S		MP03703	
Demo pack BabyFlow® plus, size M		MP03704	
Demo pack BabyFlow® plus, size L		MP03705	
Oxymixer Low Flow		MP04203	
Oxymixer Low Flow with Monitor		MP04204	
Oxymixer Low Flow with 1-15 l flowmeter right		MP04229	

Literature

¹ Short term evaluation of respiratory effort by premature infants supported with bubble nasal continuous airway pressure using Seattle-PAP and a standard bubble device. PLOS ONE, March 28, 2018, Stephen E. Welty, Craig G. Rusin, Larissa I. Stanberry, George T. Mandy, Alfred L. Gest, Jeremy M. Ford, Carl H. Backes, Jr, C. Peter Richardson, Christopher R. Howard, Thomas N. Hansen, Charles V. Smith

Not all products, features, or services are for sale in all countries.

Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53-55,
23542 Lübeck, Germany

www.draeger.com

MUMBAI (REGISTERED OFFICE)

Draeger India Private Limited
10th Floor, Commerz II,
International Business Park,
Oberoi Garden City,
Mohan Gokhale Road,
Off Western Express Highway,
Goregaon East, Mumbai - 400063
Tel +91 22 6133 5600
info.MT.india@draeger.com

BRANCH OFFICES

AHMEDABAD

2nd Floor, Panchamrut House,
Near Panchamrut - 2, Bungalows,
Near. Science City Bridge, Sola,
Ahmedabad - 380060
Tel +91 98980 00162 / 00424
Fax +91 79 2929 9101

BENGALURU

2nd Floor, No. 179, Amarjyoti Layout,
Intermediate Ring Road,
Domlur, Bengaluru - 560071
Tel +91 80 4167 2869/70
Fax +91 80 4098 5455

CHENNAI

C/o FedEx Express TSCS (I) Pvt. Ltd.
77/1, 2A, 2B, Redhills Road,
Opp. Pappad Company, Kathirvedu,
Chennai - 600066
Tel +91 44 4351 9701/02
Fax +91 44 4351 9703

KOCHI

Ground Floor, Soorya Prabha Bldg.,
Door No. 35/230D,
(40/287-A), Prabha Lane,
Automobile Road, Lane 8,
Mamangalam, Palarivattom,
Kochi - 682025
Tel +91 484 4013 118
Fax +91 484 3016 063

KOLKATA

1st Floor, Sanjoy Bhavan,
44th Shanti Pally,
Kolkata - 700107
Tel +91 33 4063 0375/76
Fax +91 33 4063 0377

LUCKNOW

101, First floor,
Shalimar Logix
Rana Pratap Marg,
Lucknow - 226001
Tel +91 0522 4006890 / 4017137
Fax +91 522 410 6890

NEW DELHI

B2-3, Ground Floor,
Africa Avenue Road,
Safdarjung Enclave,
New Delhi - 110029
Tel +91 11 4604 8436/40
Fax +91 11 4615 1401

Locate your Regional Sales
Representative at:
www.draeger.com/contact

